

AqT® Biotin Goat Anti-Rabbit IgG, Fcγ Fragment Specific, Cross-Adsorbed

100 µg of conjugate per tube

Product Number: **CM86176**

Product Description

AqT® Biotin Goat Anti-Rabbit IgG is a biotinylated polyclonal secondary antibody specific to the Fc region of Rabbit IgG heavy chains. The antibody is cross-adsorbed to minimize cross-reactivity with human serum proteins.

The product is site-specifically biotinylated at reduced Fc-region disulfide bonds using CellMosaic's proprietary super-hydrophilic AqT® linker. By positioning biotin away from the antigen-binding regions, this controlled labeling approach is designed to help preserve antibody-binding activity while providing consistent and accessible biotin placement.

Each antibody contains an average of 3–5 biotin molecules, supporting strong streptavidin-mediated signal amplification and sensitive, reproducible detection across a variety of immunoassay applications.

For more information about AqT® Biotinylated Secondary Antibodies, visit:

<https://cellmosaic.com/biotinylated-secondary-antibodies/>

The product is sold as 1 tube of 100 µg of conjugate. For bulk quantities and custom packaging requirements, please contact us for pricing and availability.

Product Information


- **Chemical Name:** AqT® Biotin Goat Anti-Rabbit IgG, Fcγ Fragment Specific, Cross-Adsorbed
- **Molecular Weight:** Approximately 162 kDa
- **Physical Appearance:** Clear solution in 1xPBS buffer with stabilizers
- **Storage Temperature:** -20 °C or below
- **Purity:** ≥95% monomer by SEC-HPLC (refer COA of each lot for actual value)
- **Average Biotin Loading:** 3–5 biotin molecules per antibody (refer COA of each lot for actual value)

Antibody Information

- **Target:** Rabbit IgG
- **Immunogen:** Rabbit IgG, Fcγ Fragment Specific
- **Species Reactivity:** Rabbit
- **Cross-Adsorption:** Human serum proteins
- **Host Species:** Goat
- **Class/Isotype:** Goat IgG
- **Antibody Type:** Polyclonal
- **Country of Manufacture:** USA

Key Features

- Site-specific biotinylation with an average of 3–5 biotin molecules per antibody
- ≥95% monomeric purity by SEC-HPLC

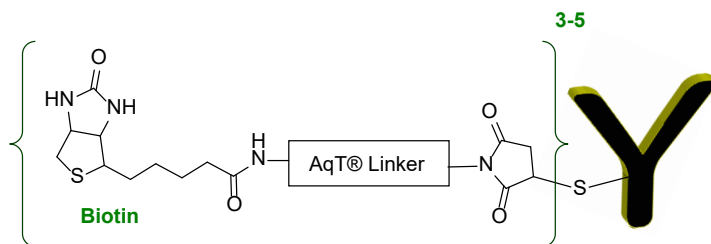

- Protein content determined by UV/HPLC analysis
- Biotin incorporation measured using a biotin assay
- Cross-adsorbed to minimize cross-reactivity with human serum proteins
- **Advanced AqT® linker technology:** CellMosaic's proprietary, super-hydrophilic, water-soluble, and charge-neutral AqT® linker is designed to support aqueous solubility and minimize aggregation


Applications

- Enzyme-linked immunosorbent assays (ELISA)
- Flow Cytometry, Immunohistochemistry (IHC)
- Immunocytochemistry (ICC)
- Western Blotting
- Immunoprecipitation (IP) and Pull-Down Assays
- Biosensor and Surface Immobilization, including SPR and BLI
- Cell Isolation and Magnetic Separation

Storage/Stability

Store the conjugate solution at -20 °C or below. Avoid repeated freeze-thaw cycles.

Figure 1. Representative size-exclusion HPLC profiles of Goat Anti-Rabbit IgG, Fcγ fragment specific before labeling (blue trace; ~2% aggregation, indicated by the peak at 10.318 minutes) and after AqT® biotin labeling (red trace). The biotinylated antibody (Lot No. S720.S9.0818F; DOL: 4.9) demonstrated ≥99% monomeric purity in PBS, with no detectable aggregate peak.

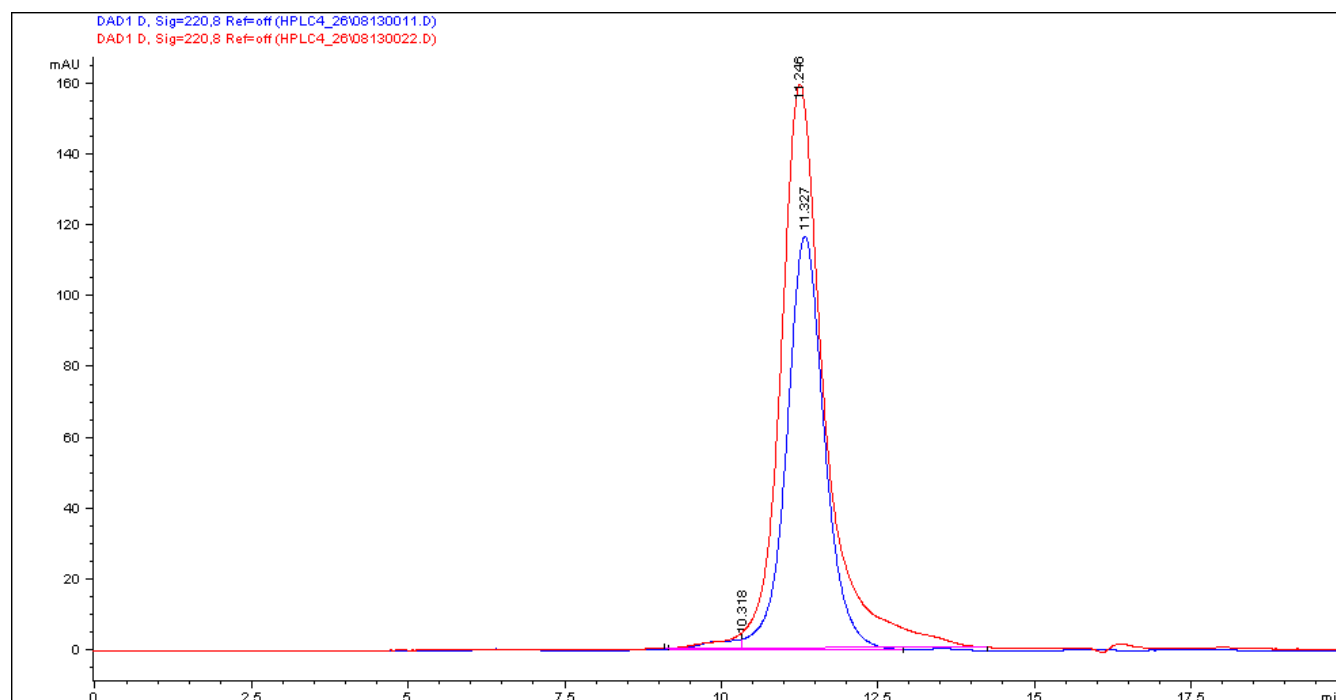

Important Notes

For Research Use Only. Not for Use in Human or Diagnostic Procedures.

The information contained in this document is subject to change without notice. CellMosaic assumes no responsibility for any errors or omissions in this document. Under no circumstances shall CellMosaic be liable, whether under contract, tort, warranty, statute, or any

other legal theory, for any special, incidental, indirect, punitive, multiple, or consequential damages related to or arising from this document or the use of this product.

Intellectual Property

CellMosaic's products, including but not limited to AqT[®], Nelon[™], oxLink[™] and sxLink[™] (collectively, the "Products"), are protected by issued patents and pending patent applications in the United States and other jurisdictions. A current list of applicable patents is available upon request.

Restricted Use

The purchase of the Product conveys only a limited, non-exclusive, non-transferable, and non-sublicensable license to use the Product for the purchaser's internal research and development purposes. No right or license is granted, whether expressly or by implication, estoppel, or otherwise, under any patent, trademark, copyright, or other intellectual property right of CellMosaic, except as expressly stated herein.

Any commercial use of the Product is strictly prohibited without a separate written agreement with CellMosaic. Prohibited commercial uses include, without limitation: (a) the sale, lease, licensing, distribution, or other transfer of the Product or any materials derived from or produced using the Product; (b) the sale, lease, licensing, or other grant of rights to use the Product or any materials derived from or produced using the Product; and (c) the use of the Product to perform services for third parties for a fee or other consideration, including, without limitation, contract research, screening services, and diagnostic or therapeutic applications

For information regarding commercial licensing or collaboration, please contact us at info@cellmosaic.com.